

Specifications

Model	Pro-care Auto Bariatric	
Pump	Dimensions	34 x 13.5 x 20.5 cm
	Weight	3.8 Kg
	Case Material	Flame retardant ABS
Supply Voltage	AC 220 - 240V , 50Hz	
	Operating Cycle	10, 15, 20, 25 mins
	Mattress	Bariatric 8" Replacement
Dimensions	200 x 107/122 x 20 cm	
	Cell Height	20 x 8" cells (Cell - on - Cell)
	Weight	13 kg
Cover Material	4-way stretch PU / Polyester cover, Nylon PU quilted cover (Optional)	
Cell Material	Nylon PU	
Maximum Patient Weight	450 kg	
Flame Retardant Standards	EN 597-1/597-2	

Specifications are subject to change for improvement without notice



HEADQUARTERS

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@apexmedicalcorp.com

PAN-ASIA

TAIWAN
Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5600 /0800 054 133
F +886 2 2268 6526
customer@apexmedicalcorp.com

STURDY INDUSTRIAL
168, Sec.1, Zhongxing Road, Wugu District,
New Taipei City 24872, Taiwan
T +886 2 2981 5245
F +886 2 2984 8208
foreign@sturdy.com.tw
www.sturdy.com.tw

CHINA
Apex Medical (Kunshan) Corp.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@apexmedicalcorp.com

THAILAND
Apex Medical (Thailand)
No.111/152 Soi Pho Kao Yak 19,
Nawamin, Beung Kum,
Bangkok 10240,
Thailand.
T +662 018 2140 Ext 4210
sales.thai@apexmedicalcorp.com

MIDDLE EAST • AFRICA
29, Misr Littameer Building,
Stree#10, AREA#8, Sheraton Building
Heliopolis Cairo, Egypt.
T +202 2 2671303
F +202 2 2678649
raafat.wahby@apexmedicalcorp.com

PAN-EUROPE

SPAIN
Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao, Vizcaya, Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedical.es

UK
Apex Medical Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester, WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales@apexmedicalcorp.co.uk

FRANCE
Apex Medical France
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d'Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@apexmedicalcorp.com

GERMANY
SLK Vertriebsgesellschaft mbH
Am Hardicksbach 18, 45731 Waltrop
Germany
T +49 (0) 231-925360-0
F +49 (0) 231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA
Apex Medical USA
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@apexmedicalcorp.com



Pro-care Auto Bariatric
Alternating Pressure Redistribution System

- Suitable for patients in high to very high risk of pressure injuries
- Total dynamic in Cell-on-Cell air cells for enhanced comfort and stable support
- Alternating Low Pressure mode for an even higher comfort level to patients
- Deliberately crafted dual compressors system silently and continuously care for patients



Distributed by

Pro-care Auto Bariatric

Alternating Pressure Redistribution System

"Pro-care, perfectly comfortable"

Built for Professionals

Pro-care Auto Bariatric is designed to meet the clinical needs in healthcare facilities to aid in pressure care plan while improving patient's comfort level, this simple and durable system helps medical professionals to remain to be highly efficient even during their challenging daily tasks.



80% of Time Saving
During decontamination, air cells with quick release detachable base can be easily removed from mattress



Max Patient Weight
Newly designed mattress system can offer high comfort and stability to patient weighs up to 450 kg

Micro Low Air Loss
Provide good ventilation and help reduce the accumulation of heat and moisture

Soft & Comfort Air Cells
A perfect selection of materials to increase comfort level for patient and reduce the noise when patient's moving on the mattress

Heel Relief Function
Heel relief function at the last 5 air cells to prevent pressure injury from occurring on patient's heels

Quiet Pump
Lower vibration from pump to enhance the better quality of patient care (<20 dBA)

A longer Lasting System Thanks to Dual Compressors
Having the capability of not only shortening mattress preparation time but also extending period of use

Graphic Design on The Panel
Easy-to-read display for a faster set-up



Non-Slip PU Base
Prevents mattress from slipping out the bed

4-way Stretch Cover
Reduced shear and low friction that further aids in pressure injury prevention. It's with welded seams that can prevent liquid from ingress into the mattress

- Low friction and shear
- Highly vapour permeable
- Water resistant
- High-flame retardant



4-Way Stretch Vapour Permeable Flame Retardant Water Resistance



Waterproof Gadget
Preventing liquid from infiltrating inside of the mattress during decontamination





Quick Connector
Once with cap on, instantly switch to static air mattress for transportation need (12 hours)



Auto Adjustment
Automatically provides an optimal pressure based on the patient's weight and position



Tuning Function
To fine-tune the pressure & optimize patient's comfort

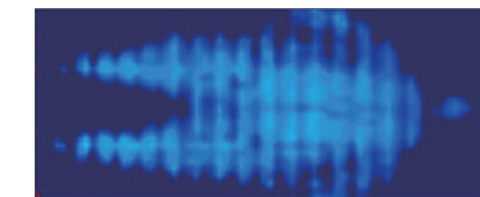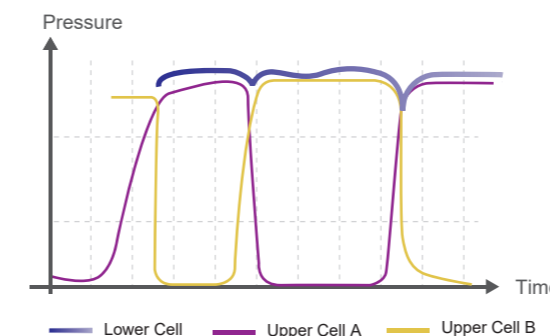

Pressure mapping under Alternating Low Pressure mode

Total Dynamic


Dynamic pressure change also in the bottom of Cell-on-Cell air cells while being able to prevent patient from bottoming out during power failure for 8 hours

Alternating Low Pressure, A Mild Alternating

Newly designed function delivers moderate alternation for sensitive patients who seek a **greater comfort experience**

Four Cycle Time Options

Alternating and Alternating Low Pressure therapy can be set at 10, 15, 20 or 25 minutes cycle based on caregiver's clinical assessment and care plan

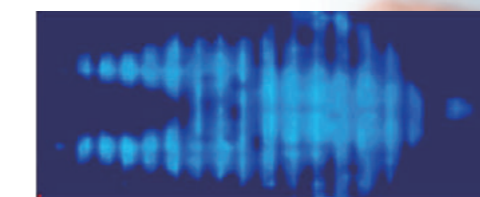

Pressure mapping under Alternating mode

